**Oral-B®**

**PRO**








powered by **BRAUN**

**TRIZONE**

**400-500**

Charger Type 3757  
Handle Type 4729

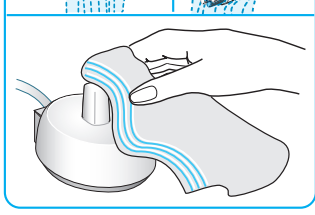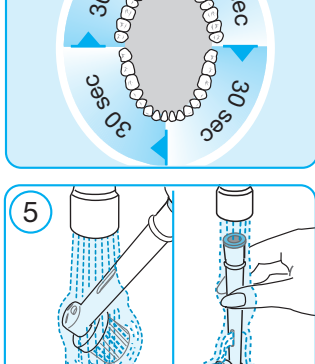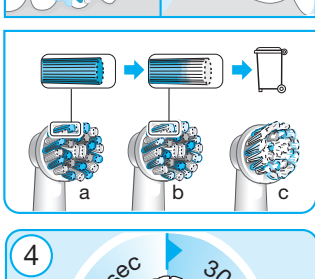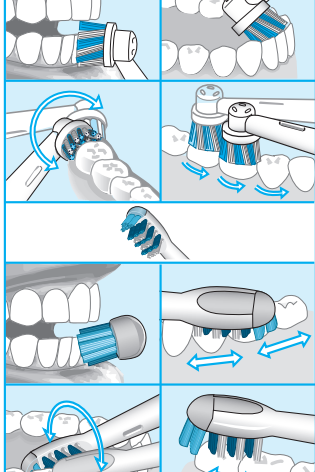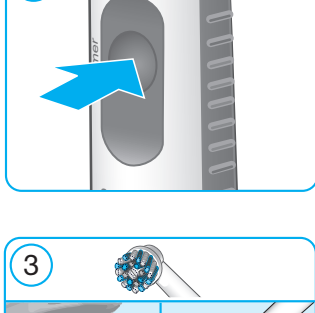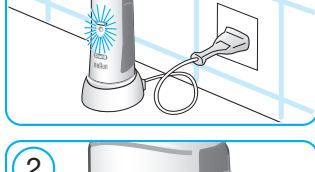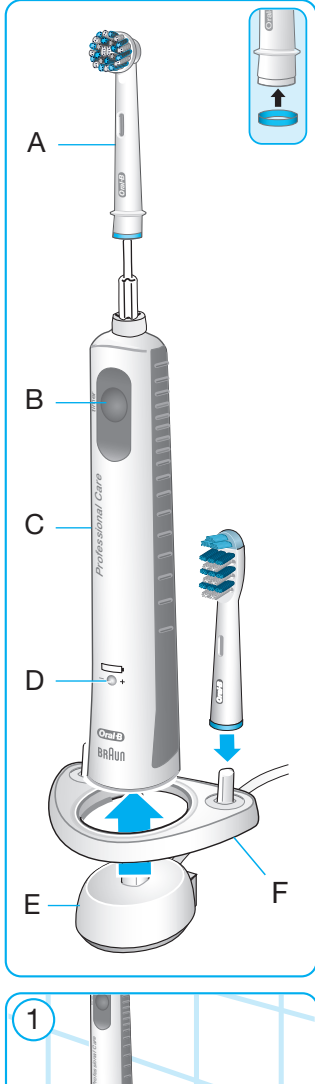- (GB)** 0 800 731 1792
- (IE)** 1 800 509 448
- (PL)** 801 127 286  
801 1 BRAUN
- (CZ)** 221 804 335
- (SK)** 02/5710 1135
- (HU)** (06-1) 451-1256
- (IL)** 1-800-666-775

www.oralb.com  
 www.braun.com  
 www.service.braun.com

**Charger type 3757**  
**Handle type 4729**

## English

Welcome to Oral-B.

Before operating this toothbrush, please read these instructions and save this manual for future reference.

### IMPORTANT

- Periodically check the cord for damage. If cord is damaged, take the charging unit to an Oral-B Service Centre. A damaged or non-functioning unit should no longer be used. Do not modify or repair the product. This may cause fire, electric shock or injury.
- This product is not intended for children 3 years and under. Children from 3 to 14 years and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge can use the toothbrush, if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved.
- Cleaning and maintenance shall not be performed by children.
- Children shall not play with the appliance.
- Use this product only for its intended use as described in this manual. Do not use attachments which are not recommended by the manufacturer.

### WARNING

- If the product is dropped, the brush head should be replaced before the next use even if no damage is visible.
- Do not place the charger in water or liquid or store where it can fall or be pulled into a tub or sink. Do not reach for it when fallen into water. Unplug immediately.
- Do not disassemble the product except when disposing of the battery. When taking out the battery for disposal of the unit, use caution not to short the positive (+) and negative (-) terminals.
- When unplugging, always hold the power plug instead of the cord. Do not touch the power plug with wet hands. This can cause an electric shock.
- If you are undergoing treatment for any oral care condition, consult your dental professional prior to use.
- This toothbrush is a personal care device and is not intended for use on multiple patients in a dental practice or institution.
- Replace brush (head) every 3 months or sooner if brush head becomes worn.

### Description

- A Brush head
- B On/off button
- C Handle
- D Charge indicator
- E Basic charger
- F Brush head holder (depending on model)

### Specifications

For voltage specifications refer to the bottom of the charging unit.

Noise level: ≤65 dB (A)

### Connecting and charging

Your toothbrush has a waterproof handle, is electrically safe and designed for use in the bathroom.

- Plug the basic charger (E) or the assembled charging unit (E/F) into an electrical outlet and place the toothbrush handle (C) on the charging unit.
- The green charge indicator light (D) flashes while the handle is being charged. Once it is fully charged, the light turns off. A full charge can take up to 22 hours and enables up to seven days of regular brushing (twice a day, 2 minutes) (picture 1). Note: After deep discharge, the green light might not flash immediately; it can take up to 10–15 minutes.
- If the rechargeable battery is low the red charge indicator light (D) is flashing for some seconds when turning your toothbrush on/off.
- For everyday use, the toothbrush handle can be stored on the charging unit to maintain it at full power. Overcharging is impossible.
- To maintain the maximum capacity of the rechargeable battery, unplug the charging unit and fully discharge the handle with regular use at least every 6 months.

### Using your toothbrush

#### Brushing technique

Wet brush head and apply any kind of toothpaste. To avoid splashing, guide the brush head to your teeth before switching on the appliance (picture 2).

When brushing your teeth with one of the Oral-B **oscillating-rotating brush heads** guide the brush head

slowly from tooth to tooth, spending a few seconds on each tooth surface.

When using the Oral-B «**TriZone**» **brush head** place the toothbrush

bristles against the teeth at a slight angle towards the gumline. Apply light

pressure and start brushing in back and forth motions, just like you would

do with a manual toothbrush.

With any brush head start brushing

the outsides, then the insides and

finally the chewing surfaces. Brush all

four quadrants of your mouth equally.

You may also consult your dentist or

dental hygienist about the right

technique for you (picture 3).

In the first days of using any electric

toothbrush, your gums may bleed

slightly. In general, bleeding should

stop after a few days. Should it persist

after 2 weeks, please consult your

dentist or dental hygienist.

#### Pressure sensor

To promote optimal brushing your

toothbrush has a pressure control

feature installed. If too much pressure

is applied, the movement of the brush

head will continue but its pulsation will

stop. In addition you will also hear a

different sound while brushing.

#### Professional Timer

A short stuttering sound at 30-second

intervals reminds you to brush equally


3

all four quadrants of your mouth (picture 4). A long stuttering sound indicates the end of the professionally recommended 2-minute brushing time. The elapsed brushing time is memorised, even when the handle is briefly switched off during brushing. When pausing longer than 30 seconds the timer resets.

## Brush heads

Our **oscillating-rotating brush heads** can be used for tooth-by-tooth precision cleaning.


### Oral-B CrossAction brush head

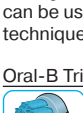Our most advanced brush head. Angled bristles for a precise clean. Lifts and powers away plaque.

Our **dynamic TriZone brush head** can be used with all common brushing techniques.

### Oral-B TriZone brush head

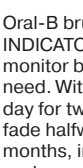features triple zone cleaning action for outstanding plaque removal, even in between teeth.

Oral-B brush heads feature light blue INDICATOR® bristles to help you monitor brush head replacement need. With thorough brushing, twice a day for two minutes, the blue color will fade halfway approximately within 3 months, indicating the need to exchange your brush head. If the bristles splay before the color recedes, you may exert too much pressure on teeth and gums (picture 3).

We do not recommend using the Oral-B «FlossAction» or Oral-B «3D White» brush head with braces. You may use the Oral-B «Ortho» brush head, specifically designed to clean around braces and wires.

## Cleaning recommendations

After brushing, rinse your brush head for several seconds under running water with the handle switched on. Switch off the handle and remove the brush head. Rinse both parts separately under running water then wipe them dry before reassembling the toothbrush. From time to time also clean the charging unit with a damp cloth; never place the charging unit in water (picture 5).

Brush head holder (F) is dishwasher-safe.

Subject to change without notice.

## Environmental Notice

Product contains batteries and /or recyclable electric waste. For environment protection do not dispose in household waste, but for recycling take to electric waste collection points provided in your country.

## Guarantee

We grant 2 years limited guarantee on the product commencing on the date of purchase. Within the guarantee period we will eliminate, free of charge, any defects in the appliance resulting from faults in materials or workmanship, either by repairing or replacing the complete appliance as we may choose. This guarantee extends to every country where this appliance is supplied by Braun or its appointed distributor.

This guarantee does not cover: damage due to improper use, normal wear or use, especially regarding the brush heads, as well as defects that have a negligible effect on the value or operation of the appliance.

The guarantee becomes void if repairs are undertaken by unauthorized persons and if original Braun parts are not used.

To obtain service within the guarantee period, hand in or send the complete appliance with your sales receipt to an authorised Oral-B Braun Service Centre.

## Warranty replacement brush heads

Oral-B warranty will be voided if electric rechargeable handle damage is found to be attributed to the use of non-Oral-B replacement brush heads. Oral-B does not recommend the use of non Oral-B replacement brush heads.

• Oral-B has no control over the quality of non-Oral-B replacement brush heads. Therefore, we cannot ensure the cleaning performance of non-Oral-B replacement brush heads, as communicated with the electric rechargeable handle at time of initial purchase.

• Oral-B cannot ensure a good fit of non-Oral-B replacement brush heads.

• Oral-B cannot predict the long-term effect of non-Oral-B replacement brush heads on handle wear.

All Oral-B replacement brush heads carry the Oral-B logo and meet Oral-B's high quality standards.

Oral-B does not sell replacement brush heads or handle parts under any other brand name.

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